


# Attends® Pull-Ons Discreet



Attends - The system with 10 absorption levels for bladder weakness

Attends Pull-Ons Discreet are a range of all-in-one absorbent pants that can be pulled up and down like normal underwear and are designed to manage very light urinary incontinence.

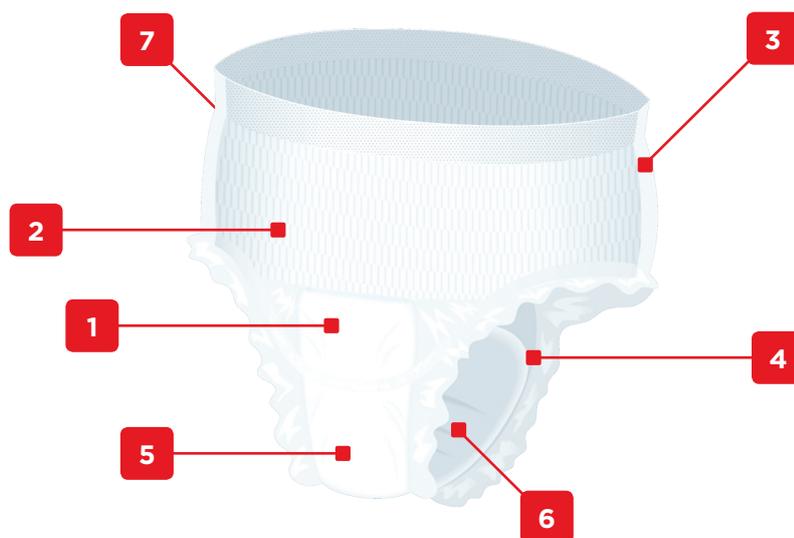

- 1** White textile feel back sheet for improved comfort and discretion.
- 2** Wide breathable elasticated waistband helps to provide a good fit so reducing the risk of leakage and improved comfort for the wearer.
- 3** Side panels which open easily if required for easy removal and disposal of the pant.
- 4** Curved leg elastics for improved fit and to help reduce the risk of leakage.
- 5** Anatomically shaped core for a close comfortable fit.
- 6** Superabsorbent core for improved leakage protection, skin dryness and odour protection.
- 7** The product size is printed on the back sheet for easy identification of products.

For more information visit [www.attends.co.uk](http://www.attends.co.uk)

**Attends Pull-Ons Discreet**

		Packcount Casecount	EAN 7332152 Pack / Case	Hip / Waist Size (cm)	Backsheet	Absorption Core Material	Topsheet Material	Elastics	Natural Odour Protection	Breathable	Quick-Dry
3	Attends Pull-Ons Discreet Medium 	12 8	205020 205037	75 - 110	Laminate of polypropylene nonwoven and liquid impervious polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓	✓
	Attends Pull-Ons Discreet Large 	10 8	205044 205051	95 - 125	Laminate of polypropylene nonwoven and liquid impervious polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓	✓



**Quick-Dry**

Quick-Dry is the acquisition layer placed directly below the top sheet where the urine enters the pad. It quickly absorbs urine away from the skin and into the pad and prevents the urine returning to the surface. It can handle several gushes of urine without affecting how the pad works which results in exceptional skin dryness and leakage protection for the user.



**Natural Odour Protection**

Natural Odour Protection reduces the risk of unpleasant smells.



**Breathable**

The breathable back sheet enables the skin to breathe, which contributes to improved skin health.