


Attends® Contours Regular



Attends - The system with 10 absorption levels for bladder weakness

Attends Contours Regular 4 - 5 are a range of anatomically shaped pads designed to manage moderate urinary incontinence.

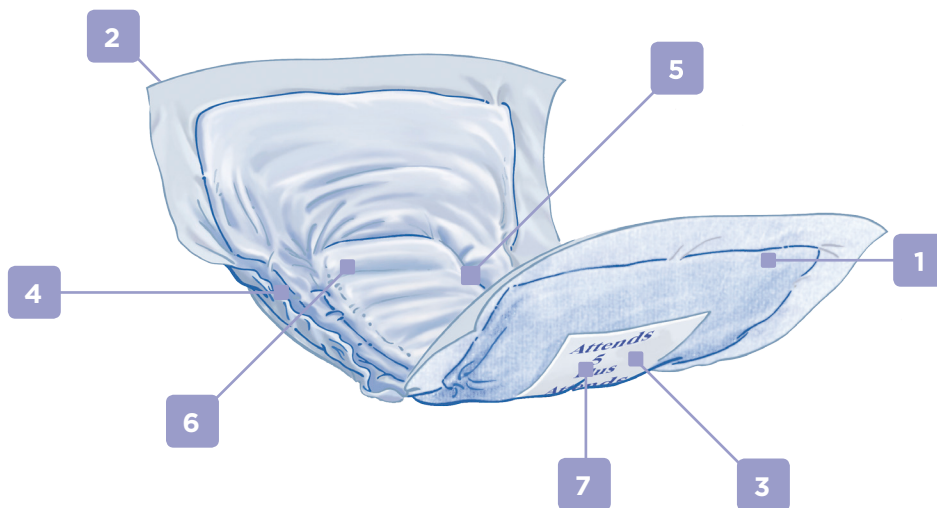

- | | | |
|--|--|--|
| 1 White waterproof back sheet to aid containment. | 2 Wetness indicator aids good nursing practice and acts as a guide to when the pad requires changing. | 3 Fixation tapes allows the pads to be worn with close fitting underwear for extra discretion. |
| 4 Soft leg elastics reduces the risk of leakage. | 5 Anatomically shaped core for a close comfortable fit. | 6 Super absorbent core with acquisition layer for improved leakage protection, skin dryness and odour protection. |
| 7 The product absorbency level is printed on the fixation tape for easy identification of products in care homes/wards. | | |

For more information visit www.attends.co.uk

	Packcount Casecount	EAN 7332152 Pack / Case	Product Length (cm) Product Width in Crotch (cm)	Backsheet	Absorption Core Material	Topsheet Material	Elastics	Natural Odour Protection	Quick-Dry Technology
Attends Contours Regular 	42 5	203491 203507	39 16	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓
Attends Contours Regular 	42 5	205211 205228	50 21	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓



Quick-Dry

The Quick-Dry technology is included in the acquisition patch directly below the top sheet in the centre of the pad. It quickly absorbs urine away from the skin and into the pad and prevents the urine returning to the surface. It can handle several gushes of urine without affecting how the pad works which results in exceptional skin dryness and leakage protection for the user.



pH Value

The core technology in the product helps to guarantee a low pH which is similar to the skin. This helps to protect the skin and ensures improved skin health.



Natural Odour Protection

Natural Odour Protection reduces the risk of unpleasant smells.



Attends Contours Regular is tested and approved by proDERM, Institute for Applied Dermatological Research.