


Attends® Contours Regular





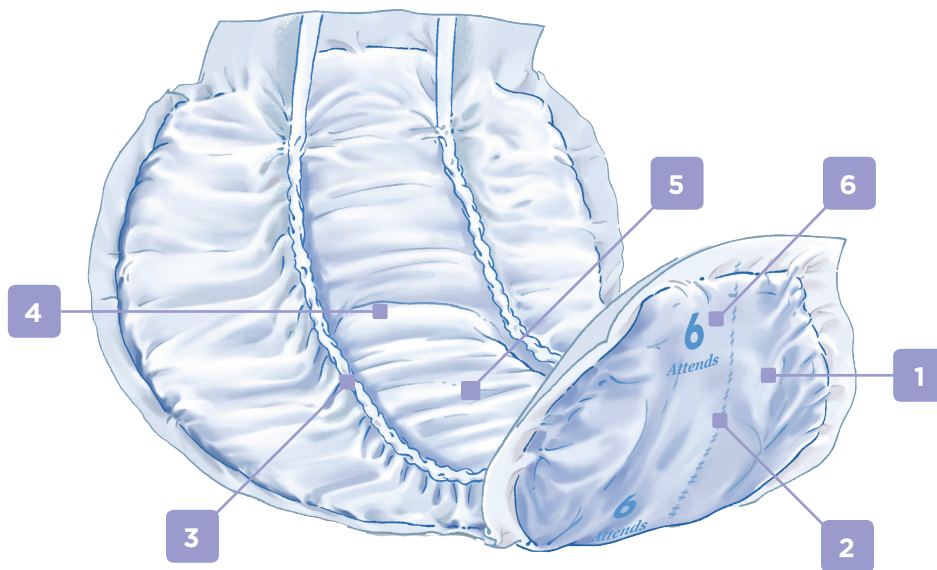


Attends - The system with 10 absorption levels for bladder weakness

Attends Contours Regular 6 - 10 are a range of anatomically shaped pads designed to manage heavy urinary and/or faecal incontinence.



- | | | |
|---|---|--|
| <p>1 White waterproof back sheet to aid containment.</p> | <p>2 Wetness indicator aids good nursing practice and acts as a guide to when the pad requires changing.</p> | <p>3 Full barrier leg cuff reduces the risk of leakage.</p> |
| <p>4 Anatomically shaped core for a close comfortable fit.</p> | <p>5 Super absorbent core with acquisition layer for improved leakage protection, skin dryness and odour protection.</p> | <p>6 The product absorbency level is printed on the back sheet for easy identification of products in care homes/wards.</p> |

For more information visit www.attends.co.uk

	Packcount Casecount	EAN 732152 Pack / Case	Product Length (cm) Product Width in Crotch (cm)	Backsheet	Absorption Core Material	Topsheet Material	Leakage Barriers	Natural Odour Protection	Quick-Dry Technology
Attends Contours Regular 6 	35 4	201244 201251	58 17.5	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓
Attends Contours Regular 7 	28 4	201268 201275	58 17.5	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓
Attends Contours Regular 8 	28 4	201282 201299	64.5 19	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓
Attends Contours Regular 9 	28 4	201305 201312	70 19.5	Polyethylene film	Fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓
Attends Contours Regular 10 	21 4	205525 205532	78.5 25	Polyethylene film	Dual Core system with fluff and superabsorbents with Quick-Dry acquisition layer	Non-woven	✓	✓	✓



Quick-Dry

The Quick-Dry technology is included in the acquisition patch directly below the top sheet in the centre of the pad. It quickly absorbs urine away from the skin and into the pad and prevents the urine returning to the surface. It can handle several gushes of urine without affecting how the pad works which results in exceptional skin dryness and leakage protection for the user.



pH Value

The core technology in the product helps to guarantee a low pH which is similar to the skin. This helps to protect the skin and ensures improved skin health.



Natural Odour Protection

Natural Odour Protection reduces the risk of unpleasant smells.



Attends Contours Regular is tested and approved by proDERM, Institute for Applied Dermatological Research.